MEDKids® 676 Baby Board, 677 Pediatric Sleeve

Model 676

Model 677


	Model 676	Model 677
		
A	23 in (584 mm)	39.75 in (1010 mm)
B	11 in (279 mm)	19.5 in (495 mm)
C	1.25 in (32 mm)	< 1 in (< 25 mm)

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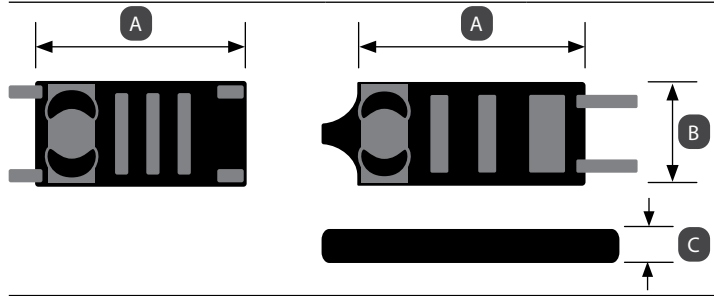Ferno (UK) Ltd.
Ferno House, Stubs Beck Lane, Cleckheaton,
West Yorkshire, BD19 4TZ, UK
+44 (0) 1274 851 999

Description

The Ferno® MEDKids® Model 676 Baby Board and Model 677 Pediatric Sleeve are patient-handling devices designed for professional use by trained operators.

Baby board is for infants 2-15 lb. (0.9-6.8 kg). It is used to restrain, position, or immobilize an infant during transportation in an isolette or on an ambulance cot.

Pediatric sleeve is for children 12-60 lb. (5.4-27.2 kg). It attaches to any backboard and is used to immobilize the patient during transportation to, and on, an ambulance cot.

General Guidelines for Use

Medical advice is beyond the parameters of this manual. It is the operators' responsibility to ensure safe practices during use of the product. A minimum of two trained operators is required. Operators must:

- Read and understand this manual.
- Train and practice on proper use of the product following standard emergency patient-handling procedures. Keep training records.
- Inspect the product and confirm it is in good condition.
- Communicate with one another and use coordinated movements to maintain control of the product.
- Stay with the patient at all times.
- Use additional trained help when working with heavy loads (patient plus equipment).
- Use as many body straps as possible to restrain the patient. Follow local protocols for positioning the patient and for strap and head block use.
- After securing the patient on the product, use the cot's restraints to secure the product on the ambulance cot before transport.
- Inflate the bladder to assist with patient positioning or airway management if permitted by local protocols.

Disinfecting and Cleaning

Remove the pediatric sleeve from the backboard before cleaning and disinfecting.

Wipe all surfaces with disinfectant. Follow the disinfectant manufacturer's instructions for application method and contact time.

Hand clean with warm water and a mild detergent. Rinse with warm, clear water. Dry with a towel or allow it to air-dry.

Ferno Customer Relations

For ordering assistance or general information, contact your Ferno distributor or Ferno Customer Relations.

Information www.ferno.com
Authorized Service ferno.com/service-and-support

U.S., Canada	1.877.733.0911
International	+1.937.382.1451
Fax	1.888.388.1349

Complete English Users' Manual
available at: www.ferno.com.

PATENTS
ferno.com/patents

WARNING

Professional use by trained operators only. Use the product only for its intended purpose.

Improper use, operation, or maintenance can cause injury. Use, operate, and maintain only as instructed.

An unattended or unrestrained patient can be injured. Always stay with the patient and use patient restraints.

Improper parts, modifications, or items not designed for the product can cause injury. Use only Ferno-approved parts and service. Attach only approved items. Do not modify the product.

Symbol	Definition	Symbol	Definition
	Follow training protocol		Load capacity
	Disinfecting		Product weight
	Cleaning		676 Baby Board
			677 Pediatric Sleeve

NOTICE

Do not use disinfectants or cleaners containing phenolics or iodines. If bleach is used, rinse with plenty of clear water to remove bleach residue.

Do not use abrasive materials to clean the product.

Storing the product wet will encourage the growth of mold. Do not store the product until it is thoroughly dry.

Inspection Checklist

- Is product free of excessive wear, punctures, frayed edges or tears?
- Are all components present?
- Is patient-restraint webbing in good condition with no cuts or frayed edges?
- Are patient-restraint buckles free of visible damage and do they operate properly?

If inspection shows damage or wear, remove the product from service and repair or replace.

Part #	Description
081-9981	Model 676 head blocks
081-9986	Model 677 head blocks
081-9983	Model 676 bladder

Adverse Event Notice

In the event of an adverse event or serious incident related to the use of this device, the end user/operator must report the incident to Ferno at 1-877-733-0911 or via email at tscoordinator@ferno.com. If the incident occurred outside of the U.S., report it to Ferno's Authorized Representative for that area.

Disposal

Prior to disposal, clean and disinfect the product according to instructions in this manual. Recycle or dispose of the product following local protocols and regulations. Do not dispose of the product in municipal waste.

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