

Winter kit

cod.21-0120-019





 **FERNO**[®]

Italy
T: +39 051 68 600 28
info.it@ferno.com
www.ferno.it

Switzerland
T: +41 41 259 60 00
info.ch@ferno.com
www.ferno.it

Winter Kit

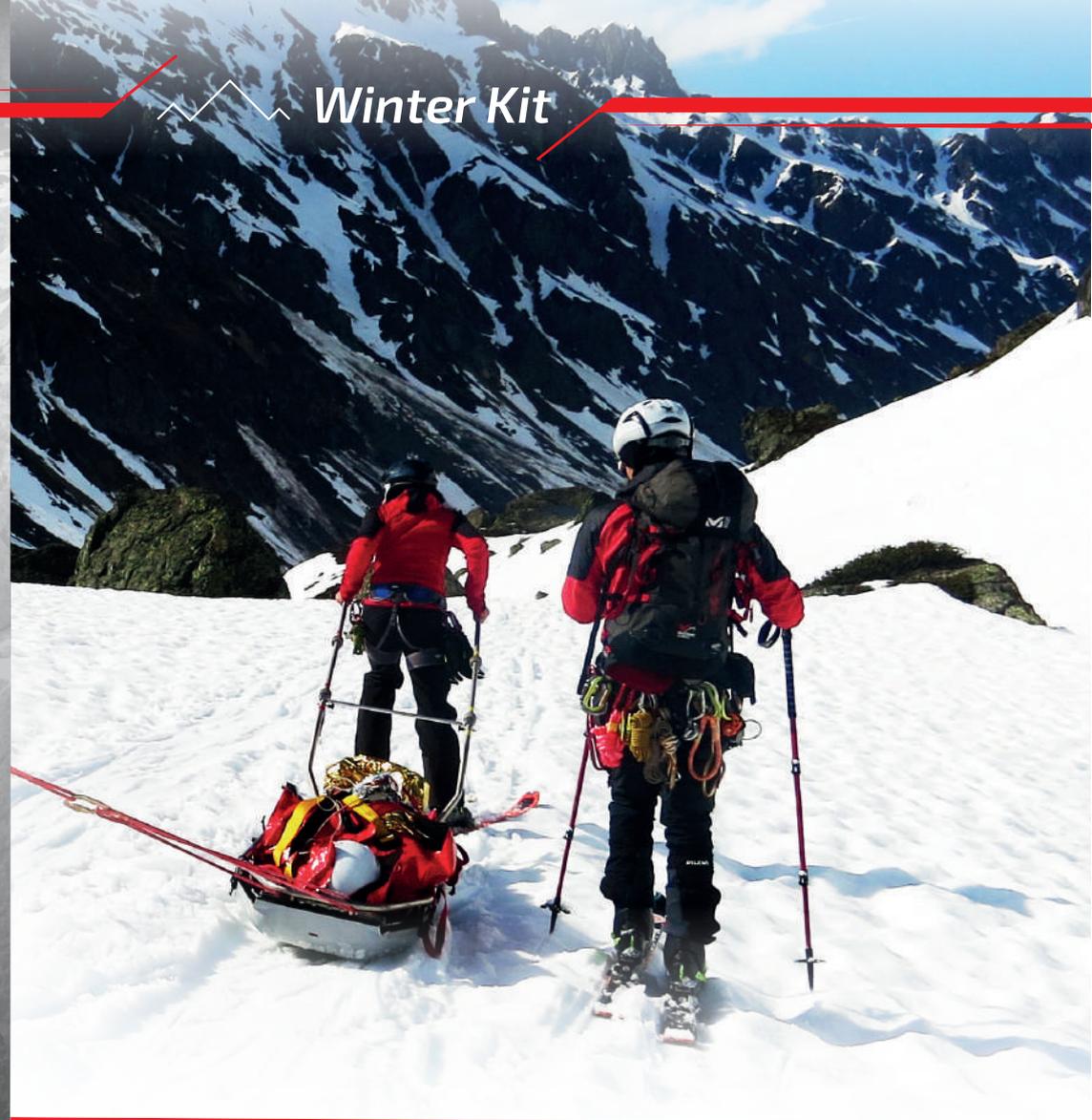

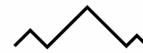Winter Kit
 For Tirol Kit

 **FERNO**[®]

Rel. 270525

This is what we call
WINTER KIT.

The rigid shells you see here make the stretcher multifunctional and suitable also in the winter.

By adding the extension handles, the stretcher can be dragged on snow and therefore you can use it on the slopes and off-piste rescue operations.

The fin-like device under the stretcher makes it stable and manageable on snow.

The shells are made of fiberglass with a kevlar inner layer to increase resistance to abrasion.

- ▶ **EASY TO INSTALL WITH SIMPLE RUBBER BANDS.**
- ▶ **COMPATIBLE WITH DIVISIBLE TITAN-T STRETCHER.**
- ▶ **THREE PIECES DIVISIBLE HANDLES AND STOWED IN SMALL PLACE.**
- ▶ **SHELLS ARE EASY TO REPAIR IF BROKEN.**

Shell



▶ **Winter kit** **cod. 21-0120-019**
no.1 head end shell, no.1 foot end shell,
handle extensions - 4pcs, no.1 connection bar

Handles

