


Model 471 Flotation Collar

A	138"	351 cm
B	4"	10 cm
	5 lb	2.27 kg
	~ 45 lb	~ 19 kg

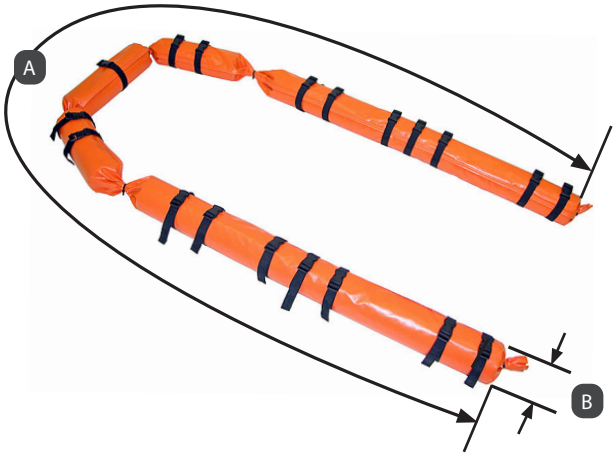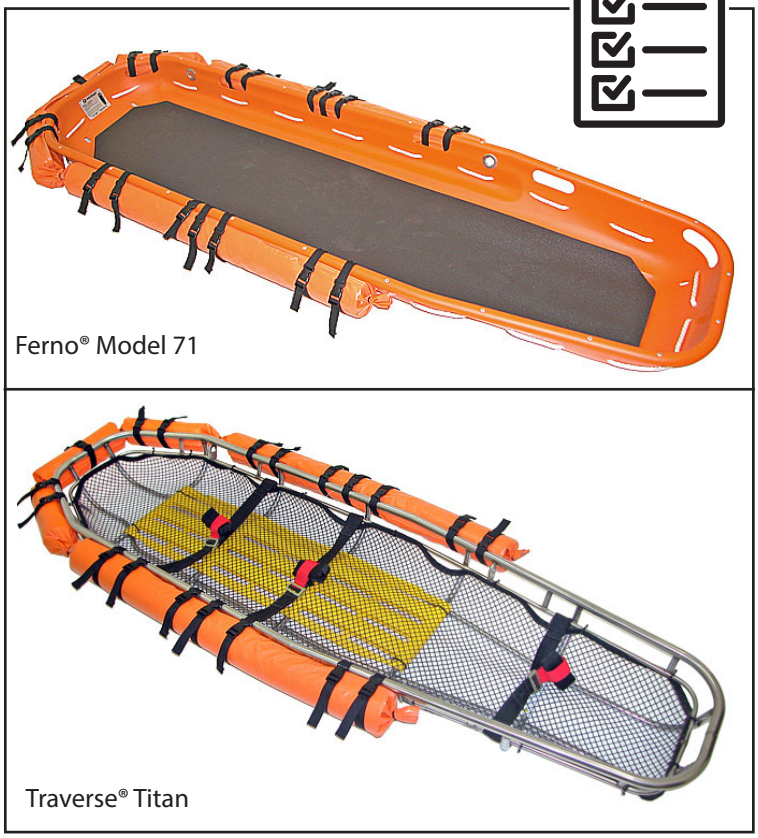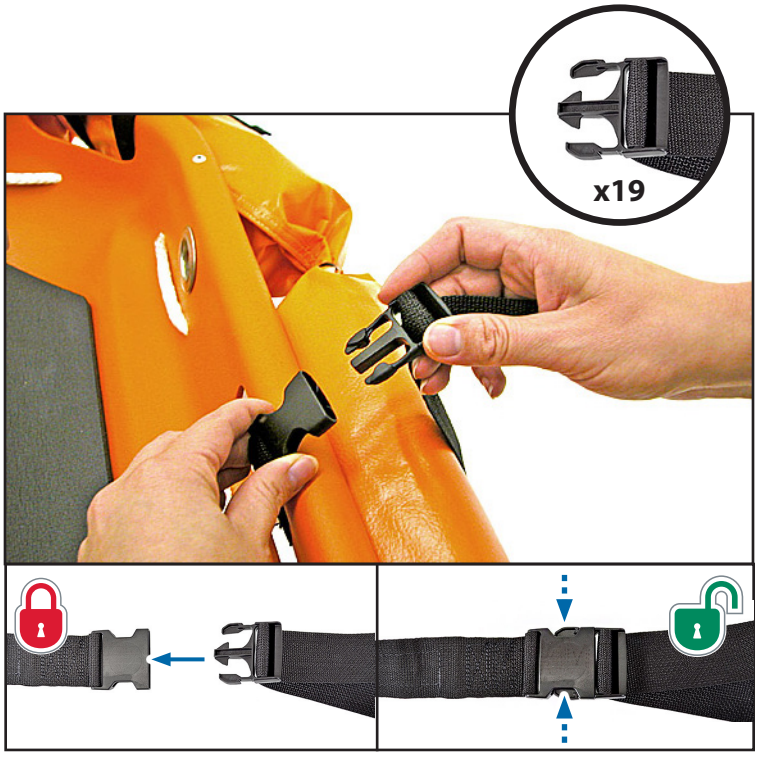Ferno-Washington, Inc.
70 Weil Way
Wilmington, OH 45177-9371, U.S.A.
1.937.382.1451

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UK REP
Ferno (UK) Ltd.
Ferno House, Stubs Beck Lane, Cleckheaton,
West Yorkshire, BD19 4TZ, UK
+44 (0) 1274 851 999

Description

The Model 471 Flotation Collar is a flotation device used to increase the buoyancy and stability of a Ferno® or Traverse® basket stretcher during a water rescue.

- Collar provides 45 lbs/19 kg of buoyancy to help keep a patient's head above water. It may be stored attached to, or inside, the basket stretcher.
- Collar is designed for use with Ferno® Model 71 series basket stretchers and Traverse® Titan™, Gazelle™, or ADVANTAGE™ series basket stretchers.
- 19 plastic-buckle fastening straps
- Polyethylene foam core
- Cover constructed of durable vinyl. The material is flame retardant and is resistant to: tearing, abrasion, ultraviolet light, shrinkage, and internally- and externally-resistant to rot and fungus (mildew).

⚠ WARNING

Professional use by trained operators only. Use the product only for its intended purpose.

Improper use, operation, or maintenance can cause injury. Use, operate, and maintain only as instructed.

An unattended or unrestrained patient can be injured. Always stay with the patient and use patient restraints.

Improper attachment can cause injury. Do not attach the collar to the patient. Attach the collar only to the basket stretcher.

Improper patient transfer can cause injury. Use two swimmers (minimum) to guide and control the basket stretcher during transfer.

Improper transport can cause injury. Do not use the collar to lift or carry the basket stretcher. Use the basket stretcher handholds for carrying and use the grommets for lifting.

Improper modification can cause injury. Do not modify the product.

NOTICE

Water rescue is a potentially hazardous activity and requires special skills. Water rescue procedures are considered technical in nature. No rescue personnel should attempt technical rescue procedures unless they have received appropriate professional instruction and training.

Personnel using the collar with a basket stretcher in a water rescue situation are responsible for knowing the techniques involved and for selecting the equipment best suited to each rescue situation.

General Guidelines for Use

Medical advice is beyond the parameters of this manual. It is the operators' responsibility to ensure safe practices during use of the product. A minimum of two trained operators is required. Operators must:

- Read and understand this manual.
- Train and practice on proper use of the product following standard emergency patient-handling procedures. Keep training records.
- Follow standard emergency water-rescue patient-handling procedures when using the collar.
- Use the collar only with a compatible basket stretcher (see list on page 1).
- Use the collar as part of a comprehensive water rescue plan.
- Use the collar only for flotation of a basket stretcher and only as described in this manual.
- Stay with the patient at all times.
- Inspect the product and confirm it is in good condition.
- Use additional trained help when working with heavy loads (patient plus equipment).

Symbol	Definition
	General warning of potential injury
	Refer to instructions for use (ANSI Z535/ISO 3864)
	Two trained operators required
	Collar weight
	Follow training protocol
	Disinfecting
	Buoyancy
	Cleaning
	Straps (19)
	Lock
	Unlock

NOTICE


Attach a towing line or lifting bridle to the stretcher main frame or around both collar and stretcher frame, never to the collar alone. The collar is not designed to bear the weight of the basket stretcher and patient.

Operation

The collar is designed to help provide buoyancy to a basket stretcher only. Do not lift or carry the stretcher by the collar.

1. To center the collar on the stretcher, start at the head end of the stretcher and attach as many of the 19 straps as possible to stretcher frame or handholds.
2. If the stretcher has a double frame, attach straps to both bars.
3. Tighten straps by pulling the loose end of each strap.
4. Lift and carry the stretcher only by hand, grasping the handholds. If the stretcher can be carried by a lifting bridle or tow lines, attach these to the stretcher, never to the collar. The collar is not designed to carry the weight of the stretcher and patient.

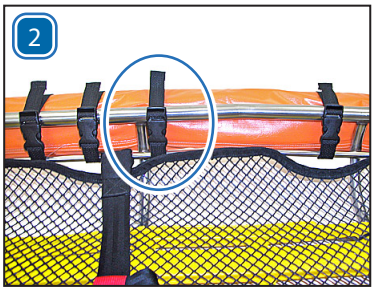Using the Basket Stretcher With the Collar

Follow your local protocols when carrying out a water rescue. Ferno recommends using a minimum of two swimmers to guide and control the basket stretcher.

1. Attach collar to stretcher before placing the stretcher in water.
2. Submerge the stretcher near the patient.
3. Guide the submerged stretcher under the patient and, while maintaining control of the stretcher, slowly bring stretcher up beneath the patient.
4. Fasten the stretcher restraints around the patient.

Disinfecting and Cleaning

Remove the collar from the stretcher before cleaning and disinfecting.

EACH USE: After each use, rinse the collar with fresh water to remove sand, salt water, etc. Do not store the collar wet. Allow the collar to air dry before storing it.

Wipe all surfaces with disinfectant. Follow the disinfectant manufacturer's instructions for application method and contact time.

Hand clean with warm water and a mild detergent. Rinse with warm, clear water. Dry with a towel or allow it to air-dry.

NOTICE

Do not use disinfectants or cleaners containing phenolics or iodines. If bleach is used, rinse with plenty of clear water to remove bleach residue.

Inspection Checklist

- Is the collar free of excessive wear, punctures, frayed edges or tears?
- Are all components present?
- Are all straps and buckles in place?
- Is strap webbing in good condition with no cuts or frayed edges?
- Are strap buckles free of visible damage and do they operate properly?

If inspection shows damage or excessive wear, remove the product from service and replace. The collar is not a repairable item.

Ferno Customer Relations

For ordering assistance or general information, contact your Ferno distributor or Ferno Customer Relations.

Information	www.ferno.com
Authorized Service and Support	ferno.com/service-and-support
U.S., Canada	1.877.733.0911
International	+1.937.382.1451
Fax	1.888.388.1349

Adverse Event Notice

In the event of an adverse event or serious incident related to the use of this device, the end user/operator must report the incident to Ferno at 1-877-733-0911 or via email at tscoordinator@ferno.com. If the incident occurred outside of the U.S., report it to Ferno's Authorized Representative for that area.

Disposal

Prior to disposal, clean and disinfect the product according to instructions in this manual. Recycle or dispose of the product following local protocols and regulations. Do not dispose of the product in municipal waste.

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© Ferno / 234-1198-03 / Febbraio 2026